

Acetaminophen capsules-USP version

Column: Luna® 5 µm C18(2) 100 Å, LC Column 250 x 4.6 mm, Ea

Dimensions: 250 x 4.6 mm ID

Order No: 00G-4252-E0

Elution Type: Isocratic

Eluent A: water


Eluent B: methanol

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
1	0	75	25	

Flow Rate: 1.5 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 243 nm (ambient)

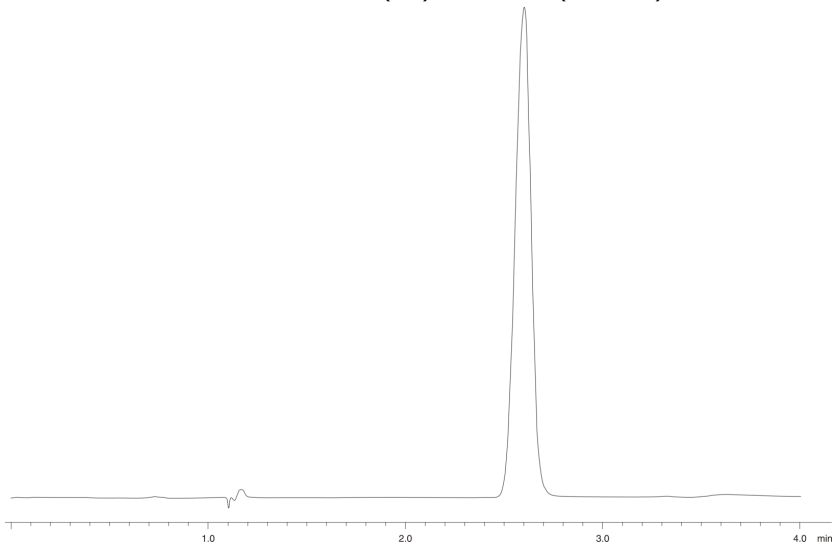

ANALYTES:

1 acetaminophen



Products used in this application:

