

Phthalates Comparison

Column: Luna® 5 µm CN 100 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4255-E0

Elution Type: Isocratic

Eluent A: Hexane

Eluent B: Methylene chloride/methanol (80:20)

Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	0	95	5

Flow Rate: 1 mL/min

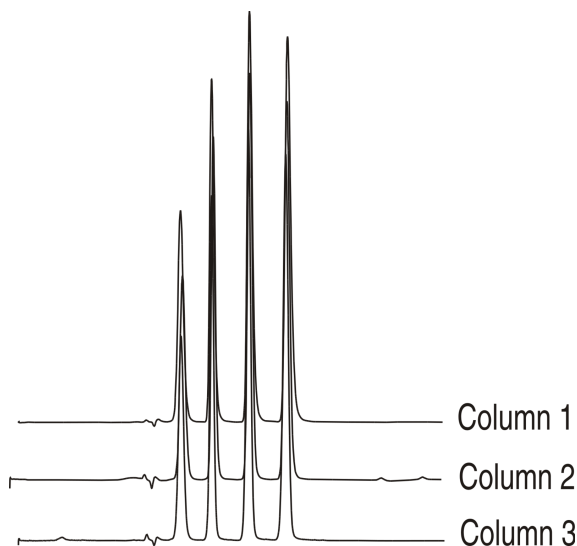
Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, CN 4 x 3.0mm, 10/Pk Part No.: AJ0-4305

- Holder Part No.: KJ0-4282

**ANALYTES:**

- 1 Dimethyl phthalate
- 2 Diethyl phthalate
- 3 Dibutyl phthalate
- 4 Dioctyl phthalate



Products used in this application:

