

Tricyclic Antidepressants

Column: Luna® 5 µm CN 100 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4255-E0

Elution Type: Isocratic

Eluent A: 20mM Potassium phosphate, pH 3.0/Methanol (80:20)

Eluent B: Acetonitrile

Gradient	Step No.	Time (min)	Pct A	Pct B
----------	----------	------------	-------	-------

Profile:	1	0	80	20
----------	---	---	----	----

Flow Rate: 1.5 mL/min


Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 210 nm (ambient)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, CN 4 x 3.0mm, 10/Pk Part No.: AJ0-4305

- Holder Part No.: KJ0-4282

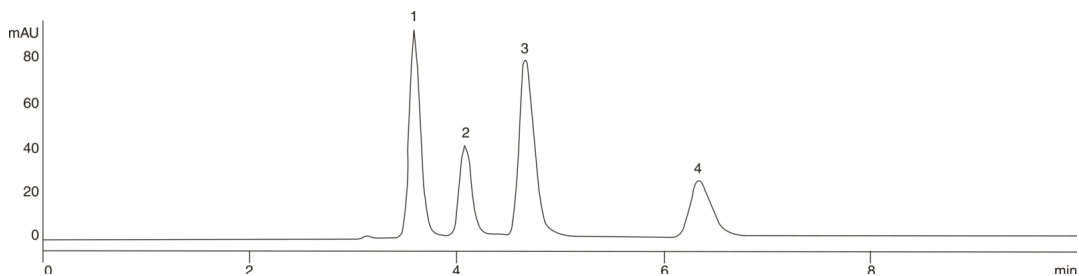

Products used in this application:

**ANALYTES:**

- 1 Doxepin
- 2 Desipramine
- 3 Amitriptyline
- 4 Clomipramine

