


EPA Method 624 Calibration Mix**Column:** Zebron[™] ZB-624, GC Cap. Column 60 m x 0.25 mm x 1.40 μ m, Ea**Phase:** 6% Cyanopropylphenyl 94% Dimethylpolysiloxane**Dimensions:** 60 meters x 0.25 mm x 1.4 μ m**Order No:** 7KG-G005-27**Oven Profile:** 40°C (hold 4 min) to 160°C @ 5°C/min**Carrier Gas:** Constant Flow Hydrogen, 1.03 mL/min**Injection:** Split 100:1 1 μ L @ 250°C**Detection:** Flame Ionization (FID) (255°C)**Products used in this application:**

EPA Method 624 Calibration Mix

ANALYTES:

- 1 1,1-Dichloroethene
- 2 Methylene chloride
- 3 1,1-Dichloroethane
- 4 Carbon tetrachloride
- 5 Benzene
- 6 Trichloroethene
- 7 1,2-Dichloropropane
- 8 2-Chlorethyl vinyl ether
- 9 1,1,2-Trichloroethane
- 10 Tetrachloroethene
- 11 Dibromochloromethane
- 12 Chlorobenzene

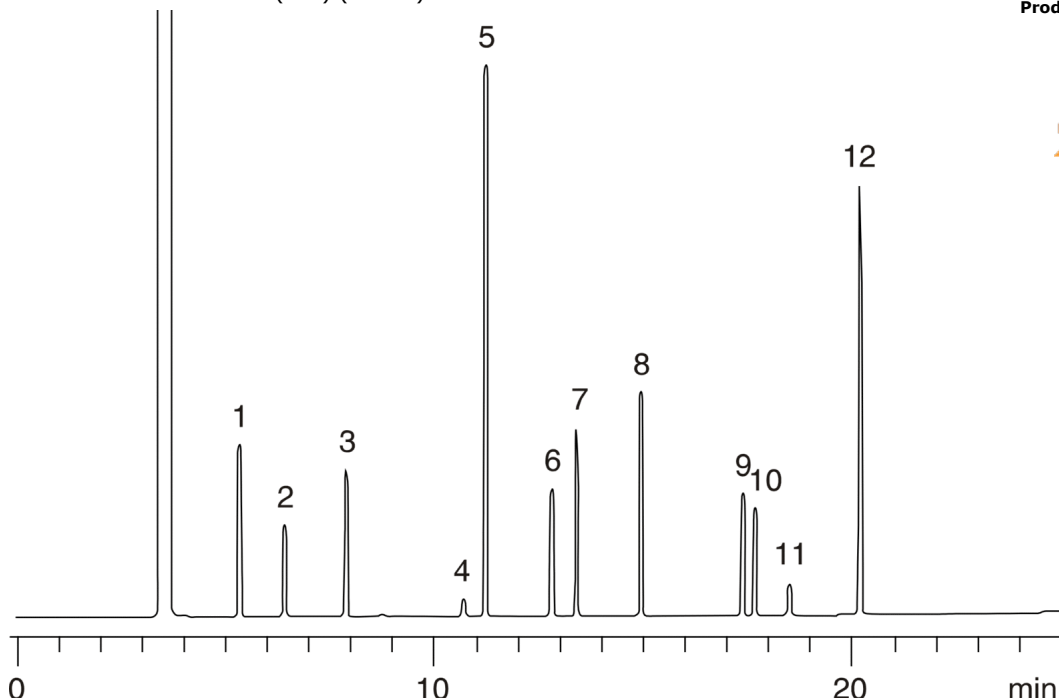