


EPA Method 608/CLP Pesticides**Column:** Zebron[™] ZB-1701, GC Cap. Column 15 m x 0.25 mm x 0.25 μ m, Ea**Phase:** 14% Cyanopropylphenyl 86% Dimethylpolysiloxane**Dimensions:** 15 meters x 0.25 mm x 0.25 μ m**Order No:** 7EG-G006-11**Oven Profile:** 100-200°C @ 20°C/min, then 5°C/min to 250°C**Carrier Gas:** Constant Flow Hydrogen, 1.47 mL/min**Injection:** Split 100:1 1 μ L @ 250°C**Detection:** Flame Ionization (FID) (285°C)**Products used in this application:**

EPA Method 608/CLP Pesticides

ANALYTES:

- 1** alpha-BHC
- 2** gamma-BHC
- 3** Heptachlor
- 4** Aldrin
- 5** beta-BHC
- 6** delta-BHC
- 7** Heptachlor epoxide
- 8** Endosulfan I
- 9** 4,4'-DDE
- 10** Dieldrin
- 11** Endrin
- 12** 4,4-DDD
- 13** Endosulfan II
- 14** 4,4'-DDT
- 15** Endrin aldehyde
- 16** Endosulfan sulfate

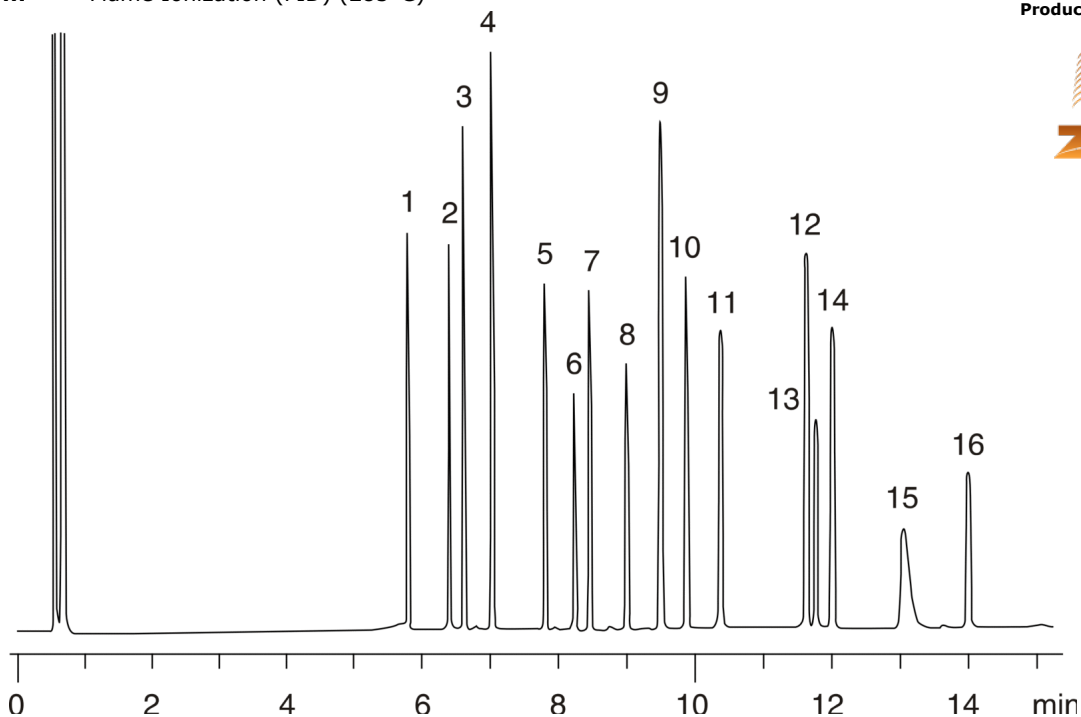