


Ylang Ylang Oil**Column:** Zebron[™] ZB-WAX, GC Cap. Column 30 m x 0.25 mm x 0.25 μ m, Ea**Phase:** Polyethylene Glycol (PEG)**Dimensions:** 30 meters x 0.25 mm x 0.25 μ m**Order No:** 7HG-G007-11**Oven Profile:** 45°C (hold 2 minutes) to 250°C @ 3°C/min**Carrier Gas:** Constant Flow Hydrogen, 1.47 mL/min**Injection:** Split 100:1 1 μ L @ 250°C**Detection:** Flame Ionization (FID) (265°C)**Products used in this application:****ANALYTES:**

- 1 Ylang Ylang oil

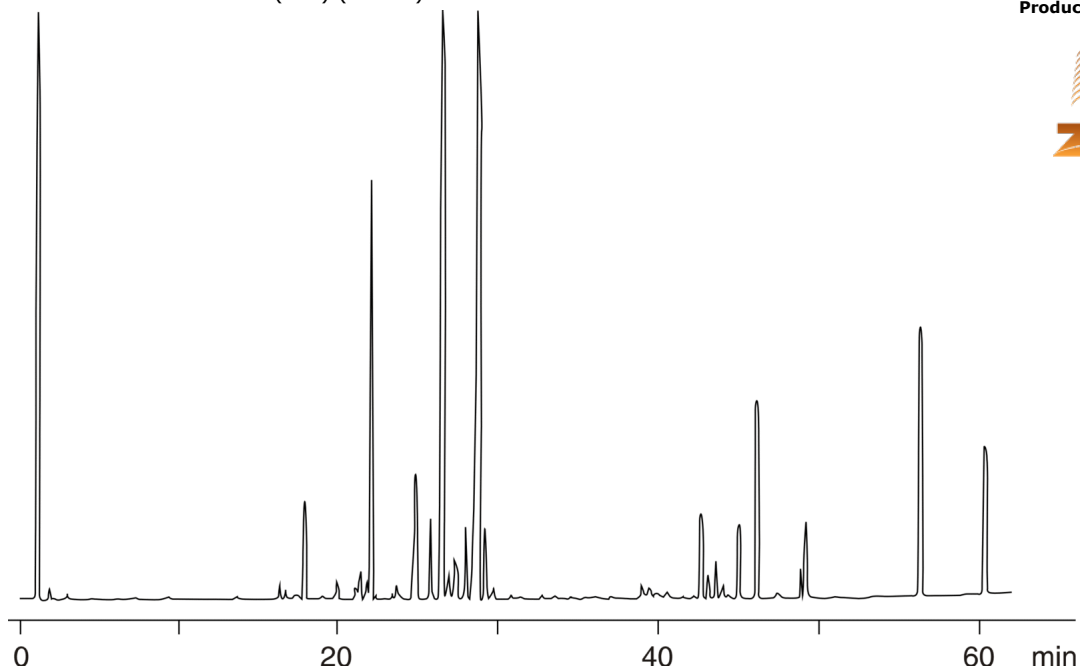