


Test Conditions for Zebron ZB-50**Column:** Zebron[™] ZB-50, GC Cap. Column 30 m x 0.25 mm x 0.25 μ m, Ea**Phase:** 50% Phenyl 50% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.25 mm x 0.25 μ m**Order No:** 7HG-G004-11**Oven Profile:** 130°C Isothermal**Carrier Gas:** Hydrogen, 1.18 mL/min**Injection:** Split 100:1 1 μ L @ 250°C**Detection:** Flame Ionization (FID) (325°C)**Products used in this application:****ANALYTES:**

- 1 Undecane
- 2 Tridecane
- 3 2,4-Dimethylphenol
- 4 2,6-Dimethylaniline
- 5 1-Undecanol
- 6 1-Methylnaphthalene
- 7 Hexadecane (C16)

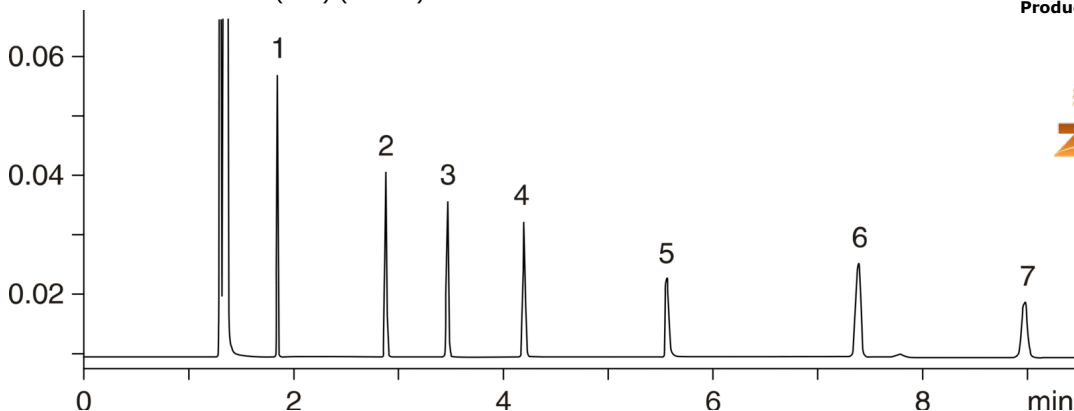