


Test Conditions for Zebron ZB-Wax**Column:** Zebron[™] ZB-WAX, GC Cap. Column 30 m x 0.25 mm x 0.25 μ m, Ea**Phase:** Polyethylene Glycol (PEG)**Dimensions:** 30 meters x 0.25 mm x 0.25 μ m**Order No:** 7HG-G007-11**Oven Profile:** 140°C Isothermal**Carrier Gas:** Hydrogen, 1.18 mL/min**Injection:** Split 100:1 1 μ L @ 250°C**Detection:** Flame Ionization (FID) (255°C)**Products used in this application:****ANALYTES:**

- 1 2-Octanone
- 2 Tetradecane
- 3 1-Octanol
- 4 Methyl decanoate
- 5 Methyl undecanoate
- 6 1-Decanol
- 7 Methyl dodecanoate
- 8 2,6-Dimethylaniline
- 9 2,4-Dimethylphenol

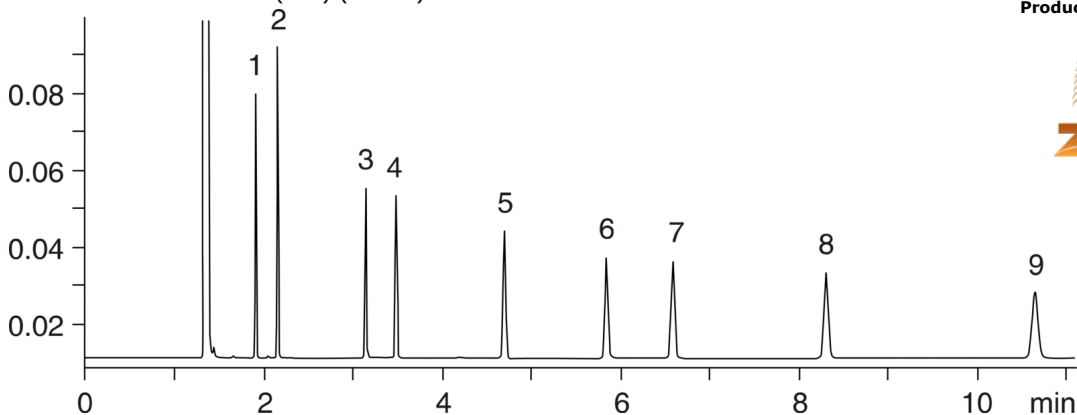