


Test Conditions for Zebron ZB-624**Column:** Zebron[™] ZB-624, GC Cap. Column 30 m x 0.25 mm x 1.40 µm, Ea**Phase:** 6% Cyanopropylphenyl 94% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.25 mm x 1.4 µm**Order No:** 7HG-G005-27**Oven Profile:** 90°C Isothermal**Carrier Gas:** Hydrogen, 1.18 mL/min**Injection:** Split 100:1 3 µL @ 250°C**Detection:** Flame Ionization (FID) (265°C)**Products used in this application:****ANALYTES:**

- 1 1,2-Dichloropropane
- 2 Octane
- 3 Tetrachloroethylene
- 4 Chlorobenzene
- 5 Nonane

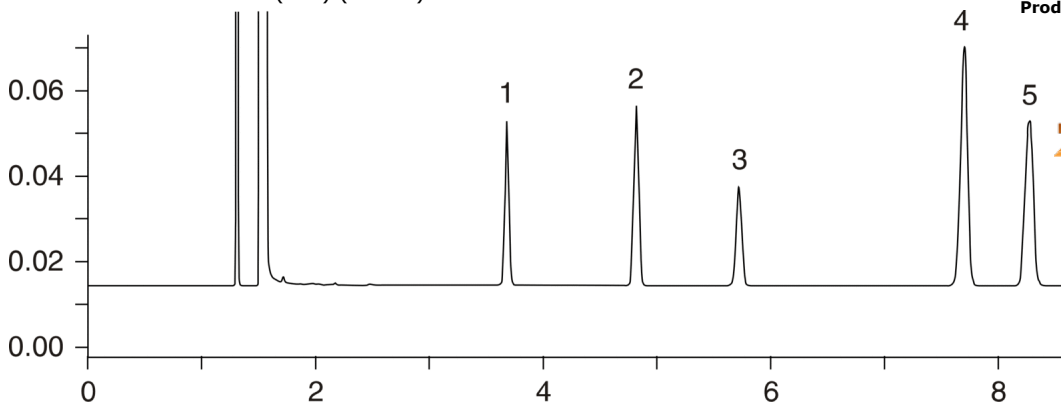