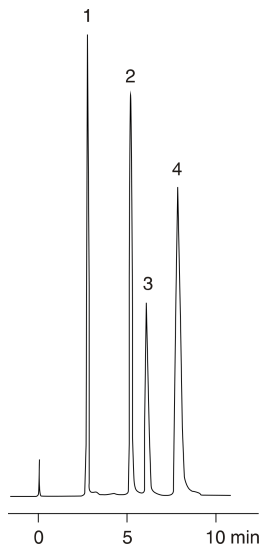


Drugs

Column: N/A
Dimensions: N/A
Order No: N/A
Elution Type: N/A
Eluent A: 25 mM Sodium phosphate pH 3.0
Eluent B: Acetonitrile
Gradient Profile:

Step No.	Time (min)	Pct A	Pct B
1	0	70	30

Flow Rate: 0.6 mL/min
Col. Temp.: ambient
Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)

**ANALYTES:**

- 1 Fluorouracil
- 2 Propranolol
- 3 Diphenhydramine
- 4 Hydrocortisone

