

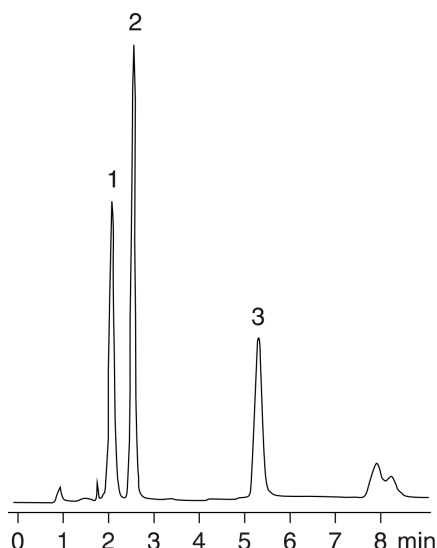
**Antioxidants - Bondclone vs.  $\mu$ Bondapak Comparison****Column:** Bondclone<sup>TM</sup> 10  $\mu$ m C18 148 Å, LC Column 300 x 3.9 mm, Ea**Dimensions:** 300 x 3.9 mm ID**Order No:** 00H-2117-C0**Elution Type:** Isocratic**Eluent A:** Acetonitrile**Eluent B:** Water

| Gradient | Step No. | Time (min) | Pct A | Pct B |
|----------|----------|------------|-------|-------|
| Profile: | 1        | 0          | 75    | 25    |

**Flow Rate:** 1.5 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)**Analyst Note:** SecurityGuard<sup>TM</sup> Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C18 4 x 3.0mm, 10/Pk Part No.: AJ0-4287

- Holder Part No.: KJ0-4282

**Products used in this application:****Bondclone<sup>TM</sup>****ANALYTES:**

- 1 TBHQ
- 2 BHA
- 3 BHT

