

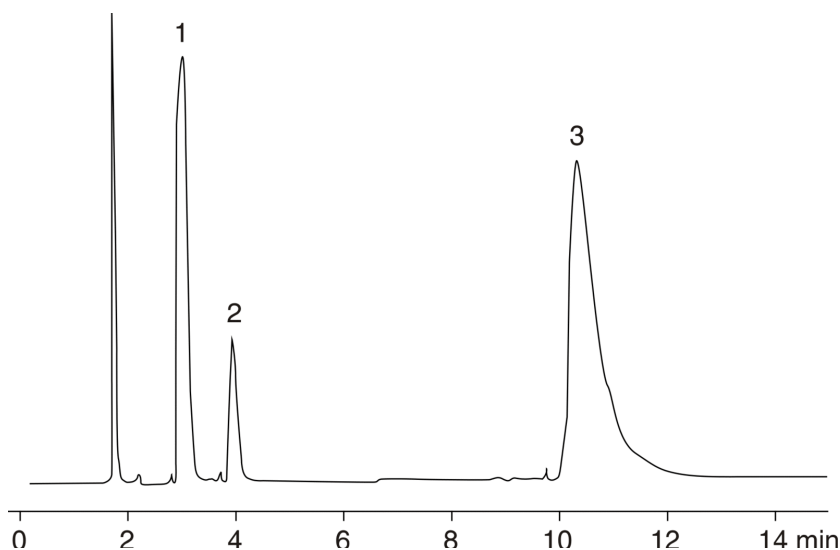
Antihistamines - Bondclone vs. μ Bondapak Comparison**Column:** Bondclone[™] 10 μ m C18 148 Å, LC Column 300 x 3.9 mm, Ea**Dimensions:** 300 x 3.9 mm ID**Order No:** 00H-2117-C0**Elution Type:** Isocratic**Eluent A:** Methanol**Eluent B:** HAC w/ 0.005 M OSA

Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	0	58	42

Flow Rate: 1.5 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)**Analyst Note:** SecurityGuard[™] Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C18 4 x 3.0mm, 10/Pk Part No.: AJ0-4287

- Holder Part No.: KJ0-4282

**Products used in this application:****Bondclone[™]****ANALYTES:**

- 1 Phenylephrine
- 2 Phenylpropanolamine
- 3 Chlorpheniramine

