

Organic Test Mix - Bondclone vs. μ Bondapak Comparison

Column: Bondclone™ 10 μ m C18 148 Å, LC Column 300 x 3.9 mm, Ea

Dimensions: 300 x 3.9 mm ID

Order No: 00H-2117-C0

Elution Type: Isocratic

Eluent A: Acetonitrile

Eluent B: Water

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	65	35

Flow Rate: 1.5 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.


- SecurityGuard Cartridges, C18 4 x 3.0mm, 10/Pk Part No.: AJ0-4287

- Holder Part No.: KJ0-4282



Products used in this application:

Bondclone™

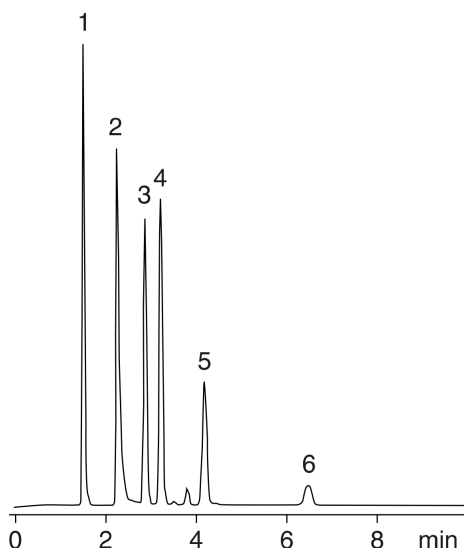

ANALYTES:

- 1 Uracil
- 2 Phenol
- 3 Benzonitrile
- 4 Nitrobenzene
- 5 1,2-Difluorobenzene
- 6 1,3-Dichloro-4,6-dinitrobenzene

