

Plant Growth Regulators - Bondclone vs. µBondapak Comparison

Column: Bondclone™ 10 µm C18 148 Å, LC Column 300 x 3.9 mm, Ea

Dimensions: 300 x 3.9 mm ID

Order No: 00H-2117-C0

Elution Type: Isocratic

Eluent A: Acetonitrile

Eluent B: Water

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	55	45

Flow Rate: 1.4 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.


- SecurityGuard Cartridges, C18 4 x 3.0mm, 10/Pk Part No.: AJ0-4287

- Holder Part No.: KJ0-4282



Products used in this application:

Bondclone™

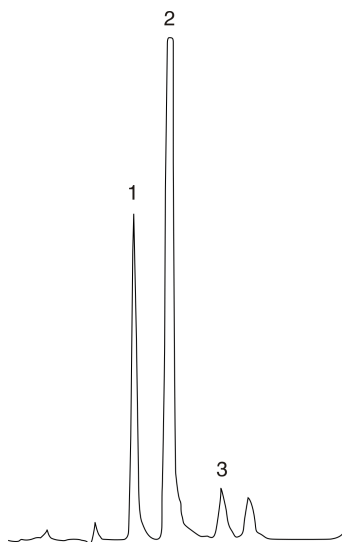

ANALYTES:

- 1 Methyl-2-chloro-9-hydroxyfluorene-9-carboxylate
- 2 Methyl-9-hydroxyfluorene-9-carboxylate

