


Methoxyverapamil**Column:** N/A**Dimensions:** N/A**Order No:** N/A**Elution Type:** N/A**Eluent A:** Hexane**Eluent B:** 1,2-Dichloroethane**Eluent C:** Ethanol-trifluoroacetic acid

Gradient	Step No.	Time (min)	Pct A	Pct B	Pct C
Profile:	1	0	55	35	10

Flow Rate: 1 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 278 nm (ambient)**ANALYTES:**

- 1 Methoxyverapamil
- 2 Methoxyverapamil

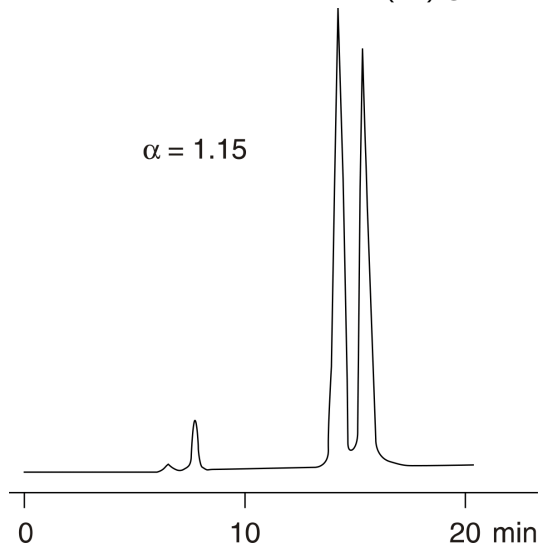