

Basic Compounds

Column: N/A

Dimensions: N/A

Order No: N/A

Elution Type: N/A

Eluent A: Acetonitrile

Eluent B: 0.01 M KH₂PO₄

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	35	65

Flow Rate: 1 mL/min


Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)

Analyst Note: SecurityGuard[™] Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C18 4 x 3.0mm, 10/Pk Part No.: AJ0-4287

- Holder Part No.: KJ0-4282

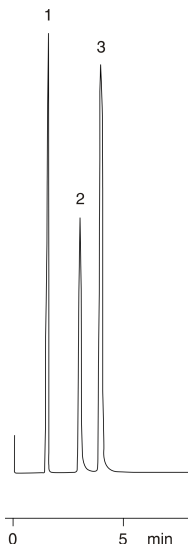

ANALYTES:

- 1 Uracil
- 2 Coptisine chloride
- 3 Berberine chloride

