


Drugs**Column:** N/A**Dimensions:** N/A**Order No:** N/A**Elution Type:** N/A**Eluent A:** Acetonitrile**Eluent B:** 0.02 M Phosphate buffer (pH 2.5)

Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	0	20	80

Flow Rate: 1 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)**Analyst Note:** SecurityGuardTM Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C18 4 x 3.0mm, 10/Pk Part No.: AJ0-4287

- Holder Part No.: KJ0-4282

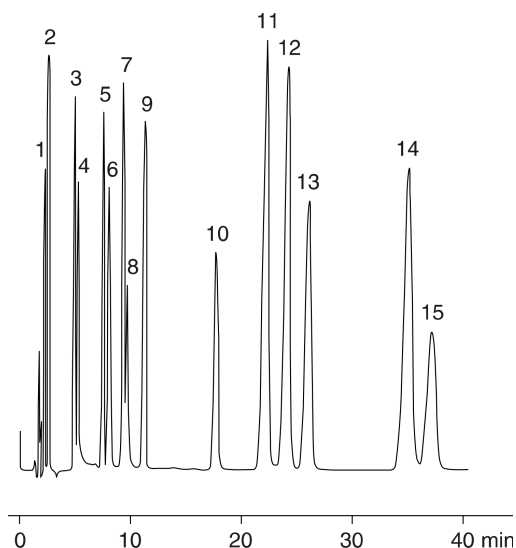

Drugs

ANALYTES:

- 1** Acetaminophen
- 2** Caffeine
- 3** Chlorpheniramine maleate
- 4** Salicylamide
- 5** Phenol
- 6** Acetylsalicylic acid
- 7** Triamcinolone
- 8** Benzoic acid
- 9** Phenacetin
- 10** Acetophenone
- 11** Prednisolone
- 12** Ethyl-p-hydroxy-benzoate
- 13** Hexobarbital
- 14** Carbamazepine
- 15** Phenytoin

