

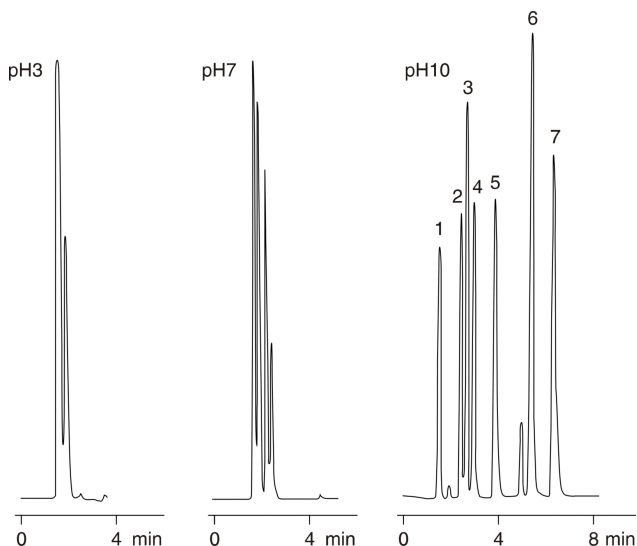
Basic Compounds**Column:** N/A**Dimensions:** N/A**Order No:** N/A**Elution Type:** N/A**Eluent A:** Acetonitrile**Eluent B:** 0.02 M Phosphate buffer

Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	0	40	60

Flow Rate: 1 mL/min**Col. Temp.:** 30 °C**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C18 4 x 3.0mm, 10/Pk Part No.: AJ0-4287

- Holder Part No.: KJ0-4282

**ANALYTES:**

- 1 Pyridoxine
- 2 Pyridine
- 3 Procainamide
- 4 Benzylamine
- 5 N-Methylbenzylamine
- 6 Metoclopramide
- 7 N,N-Dimethylbenzylamine

