

Peptide Hormones

Column: Jupiter® 5 µm C4 300 Å, LC Column 250 x 4.6 mm, Ea

Dimensions: 250 x 4.6 mm ID

Order No: 00G-4167-E0

Elution Type: Gradient

Eluent A: 0.1% TFA in water

Eluent B: 0.1% TFA in acetonitrile/water (90:10)

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	90	10
	2	8	74	26

Flow Rate: 1 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ nm (ambient)


Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Widespore C4 4 x 3.0mm, 10/Pk Part No.: AJ0-4330

- Holder Part No.: KJ0-4282



Products used in this application:

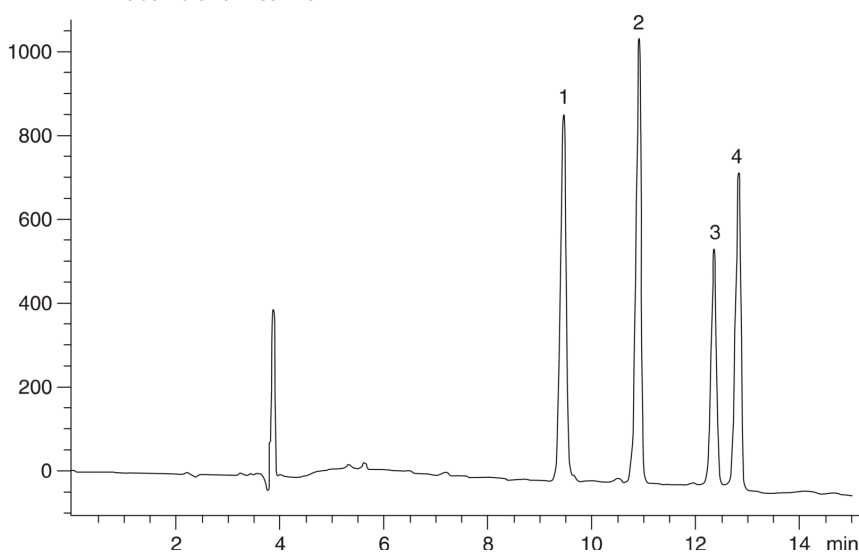

ANALYTES:

- 1 [Arg8]-Vasotocin
- 2 [Arg8]-Vasopressin
- 3 Isotocin
- 4 Oxytocin

