

Dansylated Monoamines

Column: N/A

Dimensions: N/A

Order No: N/A

Elution Type: N/A

Eluent A: Methanol

Eluent B: Water mixture

Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	0	55	45
	2	20	90	10

Flow Rate: 1 mL/min

Col. Temp.: ambient


Detection: Fluorimetric @ nm (ambient)

Detector Info: The following Fluorometric filters were employed: 7-60/7-54 and 3-72/4-76

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C8 4 x 3.0mm, 10/Pk Part No.: AJ0-4290

- Holder Part No.: KJ0-4282

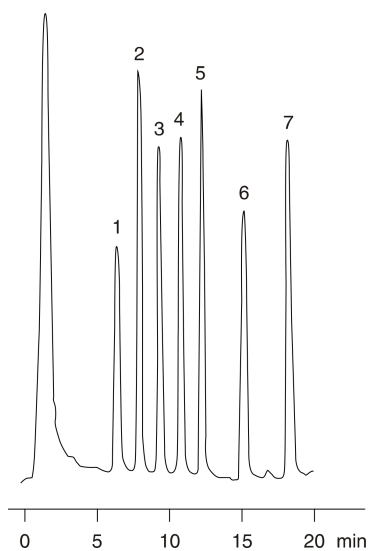

ANALYTES:

- 1 Ethylamine
- 2 Propylamine
- 3 Butylamine
- 4 Pentylamine
- 5 Hexylamine
- 6 Octylamine
- 7 Decylamine

