

Inorganic Anions

Column: N/A

Dimensions: N/A

Order No: N/A

Elution Type: N/A

Eluent A: 0.03 M Sodium salicylate, pH 4.0

Gradient **Step No.** **Time (min)** **Pct A**

Profile: **1** 0 100

Flow Rate: 0.5 mL/min


Col. Temp.: ambient

Detection: Refractive Index (RI) @ nm (ambient)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, SAX 4 x 3.0mm, 10/Pk Part No.: AJ0-4311

- Holder Part No.: KJ0-4282

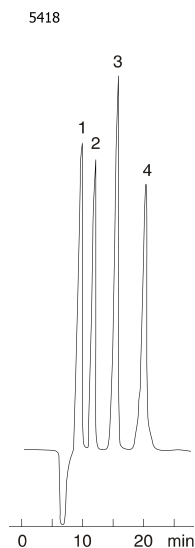

ANALYTES:

- 1** Phosphate
- 2** Chloride
- 3** Nitrate
- 4** Sulfate

