

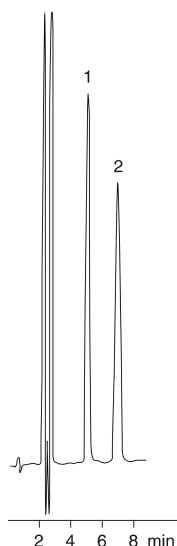
Acetanilide and Prednisone USP XXII P. 1137**Column:** N/A**Dimensions:** N/A**Order No:** N/A**Elution Type:** N/A**Eluent A:** Water**Eluent B:** Tetrahydrofuran**Eluent C:** Methanol

Gradient	Step No.	Time (min)	Pct A	Pct B	Pct C
Profile:	1	0	68.8	25	6.2

Flow Rate: 1 mL/min**Col. Temp.:** 25 °C**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ nm (ambient)**Analyst Note:** SecurityGuardTM Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, SAX 4 x 3.0mm, 10/Pk Part No.: AJ0-4311

- Holder Part No.: KJ0-4282

**ANALYTES:****1** Acetanilide**2** Prednisone