

GPC of Polystyrene Standards on Phenogel 5u 10E3Å and 10E4Å 300 x 4.6 mm columns in series

Column: Phenogel™ 5 µm 10E4A, LC Column 300 x 4.6 mm, Ea

Dimensions: 300 x 4.6 mm ID

Order No: 00H-0445-E0

Elution Type: Isocratic


Eluent A: Tetrahydrofuran

Gradient **Step No.** **Time (min)** **Pct A**
Profile: **1** **0** **100**

Flow Rate: 0.28 mL/min

Col. Temp.: 40 °C

Detection: Refractive Index (RI) @ nm (ambient)

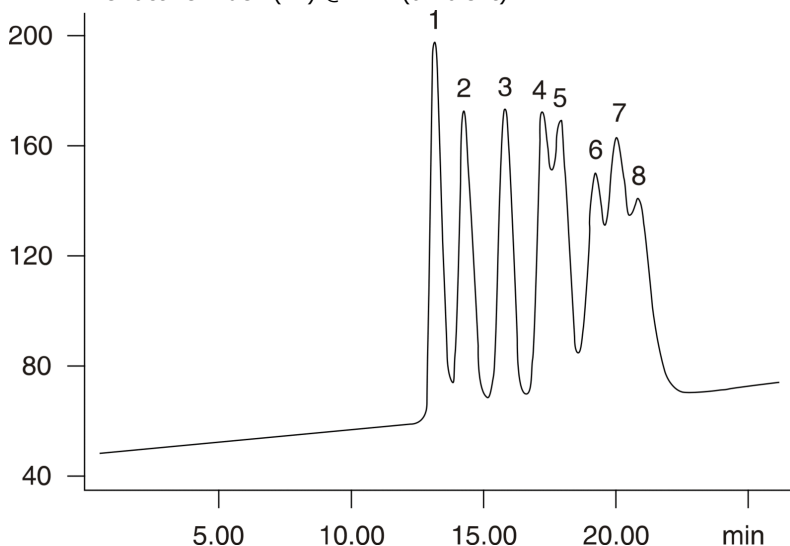

Products used in this application:



ANALYTES:

- 1** 1x10⁶ daltons
- 2** 170000
- 3** 68000
- 4** 34500
- 5** 25000
- 6** 13000
- 7** 9200
- 8** 5200

