

Polyethylene Oxide (PEO)

Column: Phenogel™ 10 µm 10E5A, LC Column 300 x 7.8 mm, Ea

Dimensions: 300 x 7.8 mm ID

Order No: 00H-0646-K0

Elution Type: Isocratic

Eluent A: DMF (0.1 M LiBr)


Gradient	Step No.	Time (min)	Pct A
Profile:	1	0	100

Flow Rate: 1 mL/min

Col. Temp.: 50 °C


Detection: Differential refractometer @ nm (ambient)

Analyst Note: 3 Phenogel 10µ columns (300 x 7.8mm ID) placed in series as follows: 10⁵Å, 10⁴Å, 10³Å.

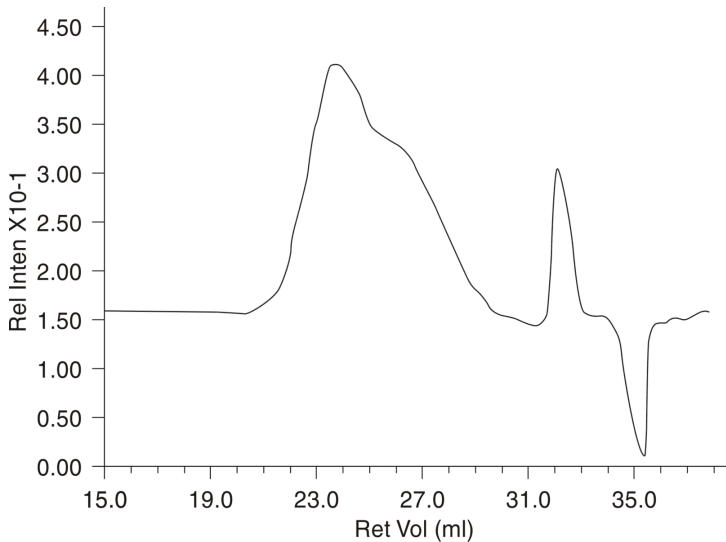

ANALYTES:

1 400,000 MW

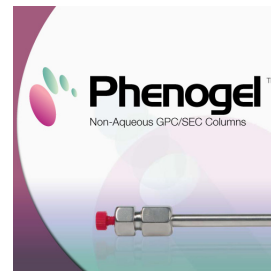

Products used in this application:

