

Mixed Polystyrene Standard on Linear(2) Column

Column: Phenogel™ 5 µm Linear(2), LC Column 300 x 7.8 mm, Ea

Dimensions: 300 x 7.8 mm ID

Order No: 00H-3259-K0

Elution Type: Isocratic


Eluent A: Tetrahydrofuran

Gradient	Step No.	Time (min)	Pct A
Profile:	1	0	100

Flow Rate: 1 mL/min

Col. Temp.: 35 °C

Detection: Refractive Index (RI) @ nm (ambient)

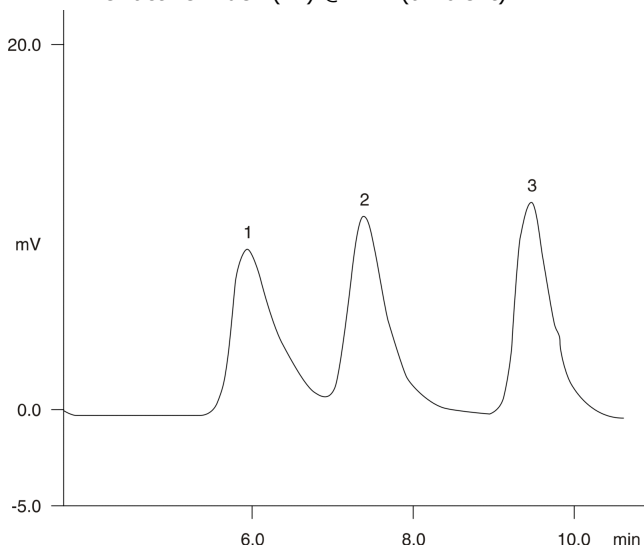

ANALYTES:

- 1 2,860,000 MW
- 2 277,000 MW
- 3 9,350 MW



Products used in this application:

