

Amitriptyline Chromatographic Comparison

Column: Prodigy™ 5 µm C8 150 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-3301-E0

Elution Type: Isocratic

Eluent A: 20 mM KH₂PO₄/KH₂PO₄, pH 7.0

Eluent B: Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	35	65

Flow Rate: 1 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C8 4 x 3.0mm, 10/Pk Part No.: AJ0-4290


- Holder Part No.: KJ0-4282

5µ C8



Products used in this application:

Prodigy™

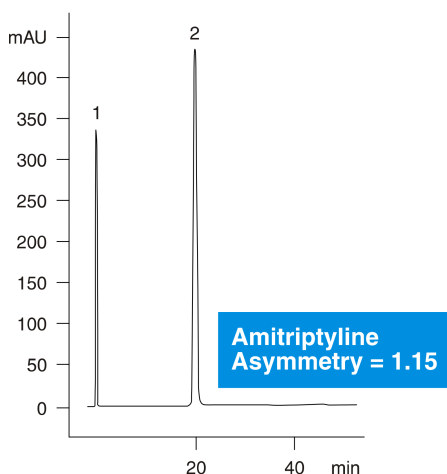

ANALYTES:

- 1 Uracil
- 2 Amitriptyline

