

Carbohydrates

Column: Rezex™ RKP-Potassium K+ 8%, LC Column 300 x 7.8 mm, Ea

Dimensions: 300 x 7.8 mm ID

Order No: 00H-3252-K0

Elution Type: Isocratic


Eluent A: Water

Gradient	Step No.	Time (min)	Pct A
Profile:	1	0	100

Flow Rate: 0.4 mL/min

Col. Temp.: 85 °C

Detection: Refractive Index (RI) @ nm (ambient)

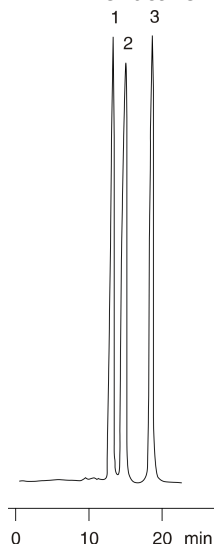

ANALYTES:

- 1 Maltotriose
- 2 Maltose
- 3 Glucose

REZEX
HPLC COLUMNS for Carbohydrates, Oligosaccharides, and Organic Acid Analysis



Products used in this application:

REZEX™

