

Oligosaccharides

Column: Rezex™ RSO-Oligosaccharide Ag+ 4%, LC Column 200 x 10 mm, Ea

Dimensions: 200 x 10 mm ID

Order No: 00P-0133-N0

Elution Type: Isocratic

Eluent A: Water

Gradient	Step No.	Time (min)	Pct A
Profile:	1	0	100

Flow Rate: 0.3 mL/min

Col. Temp.: 75 °C


Detection: Refractive Index (RI) @ nm (ambient)

REZEX™
HPLC COLUMNS for Carbohydrates, Oligosaccharides, and Organic Acid Analysis



Products used in this application:

REZEX™

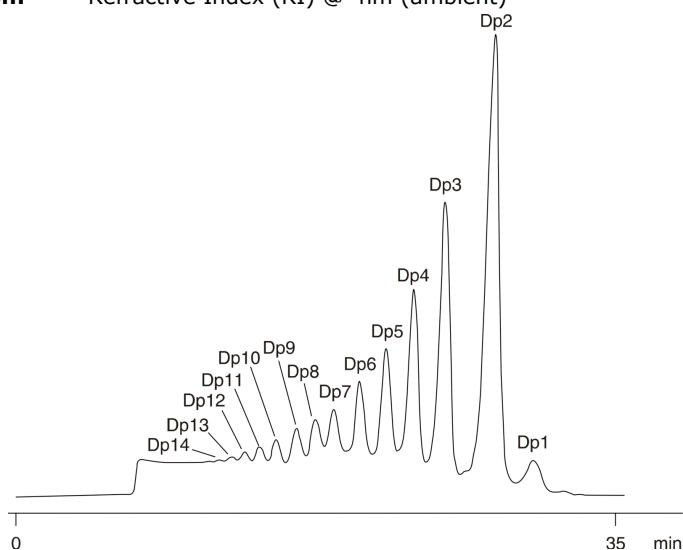

ANALYTES:

1 Dp1-Dp14

