

Carbohydrates

Column: Rezex™ RNM-Carbohydrate Na+ (8%), LC Column 300 x 7.8 mm, Ea

Dimensions: 300 x 7.8 mm ID

Order No: 00H-0136-K0

Elution Type: Isocratic


Eluent A: Water

Gradient	Step No.	Time (min)	Pct A
Profile:	1	0	100

Flow Rate: 0.6 mL/min

Col. Temp.: 75 °C

Detection: Refractive Index (RI) @ nm (ambient)

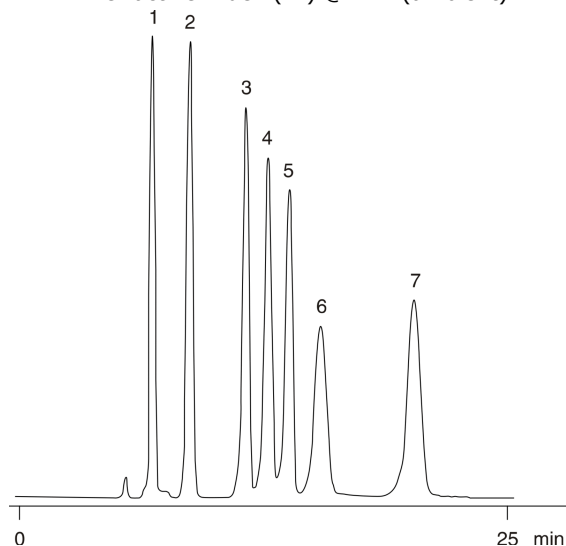

ANALYTES:

- 1 Stachyose
- 2 Cellobiose
- 3 Glucose
- 4 Fructose
- 5 Arabinose
- 6 Ribose
- 7 Salicin

REZEX™
HPLC COLUMNS for Carbohydrates, Oligosaccharides, and Organic Acid Analysis



Products used in this application:

REZEX™

