

## DMSO and Ethanol

**Column:** Rezex™ ROA-Organic Acid H+ (8%), LC Column 150 x 7.8 mm, Ea

**Dimensions:** 150 x 7.8 mm ID

**Order No:** 00F-0138-K0

**Elution Type:** Isocratic

**Eluent A:** 0.005 M H<sub>2</sub>SO<sub>4</sub>

**Eluent B:** Acetonitrile

| Gradient Profile: | Step No. | Time (min) | Pct A | Pct B |
|-------------------|----------|------------|-------|-------|
|                   | 1        | 0          | 90    | 10    |

**Flow Rate:** 0.6 mL/min


**Col. Temp.:** 40 °C

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 210 nm (ambient)

**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Carbo-H 4 x 3.0mm, 10/Pk Part No.: AJ0-4490

- Holder Part No.: KJ0-4282

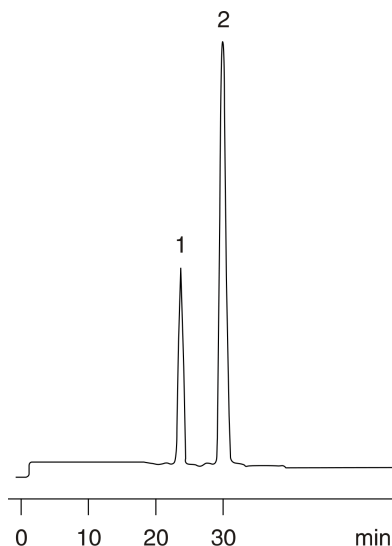

### ANALYTES:

- 1 Ethanol
- 2 Dimethyl sulfoxide

**REZEX™**  
HPLC COLUMNS for Carbohydrates, Oligosaccharides, and Organic Acid Analysis



Products used in this application:

**REZEX™**

