

4-Hydroxybenzoates - SphereClone vs. Spherisorb Comparison

Column: SphereClone™ 5 µm ODS(2) 80 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4144-E0

Elution Type: Isocratic

Eluent A: Water

Eluent B: Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
1	1	0	40	60

Flow Rate: 1 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.


- SecurityGuard Cartridges, C18 4 x 3.0mm, 10/Pk Part No.: AJ0-4287

- Holder Part No.: KJ0-4282



Products used in this application:

SphereClone™

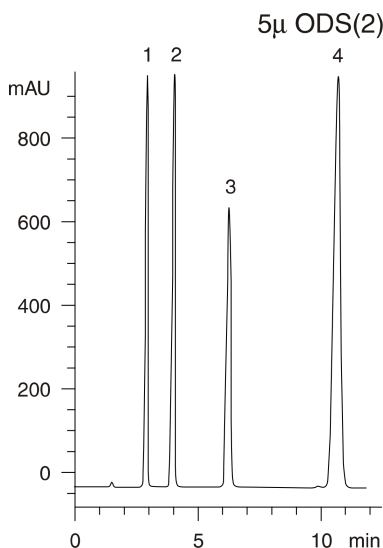

ANALYTES:

- 1 Methyl 4-hydroxybenzoate
- 2 Ethyl 4-hydroxybenzoate
- 3 Propyl 4-hydroxybenzoate
- 4 Butyl 4-hydroxybenzoate

