

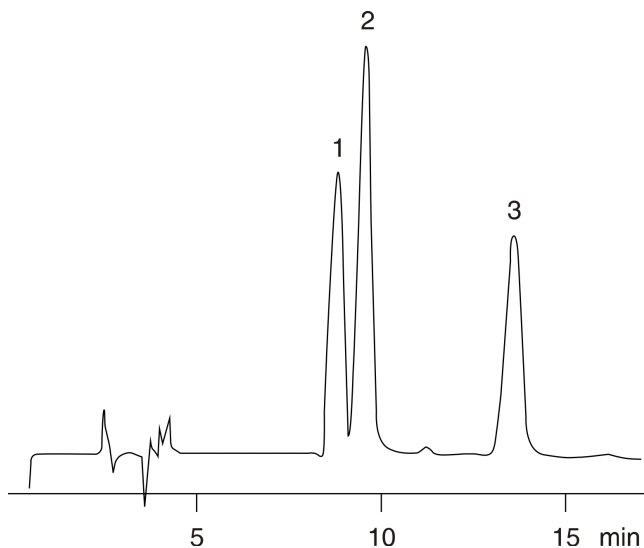
Butyrophenone Related Compound**Column:** Ultracarb[™] 5 μ m ODS (30) 60 Å, LC Column 150 x 4.6 mm, Ea**Dimensions:** 150 x 4.6 mm ID**Order No:** 00F-0351-E0**Elution Type:** Isocratic**Eluent A:** Acetonitrile**Eluent B:** 4 mM Ammonium phosphate

Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	0	70	30

Flow Rate: 0.8 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)**Analyst Note:** SecurityGuard[™] Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C18 4 x 3.0mm, 10/Pk Part No.: AJ0-4287

- Holder Part No.: KJ0-4282

**Products used in this application:****Ultracarb[™]****ANALYTES:****1** Butyrophenone