

Dabsyl Derivated Phospho Amino Acids

Column: Ultracarb™ 5 µm ODS (30) 60 Å, LC Column 250 x 4.6 mm, Ea

Dimensions: 250 x 4.6 mm ID

Order No: 00G-0351-E0

Elution Type: Gradient

Eluent A: 20 mM Sodium bicarbonate pH 8.1 with 4% dimethylformamide

Eluent B: (90:10) Acetonitrile/2-Propanol

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	82	18
	2	6	75	25
	3	26	60	40
	4	41	30	70
	5	43	30	70

Flow Rate: 1 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 460 nm (ambient)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.


- SecurityGuard Cartridges, C18 4 x 3.0mm, 10/Pk Part No.: AJ0-4287

- Holder Part No.: KJ0-4282



Products used in this application:

Ultracarb™

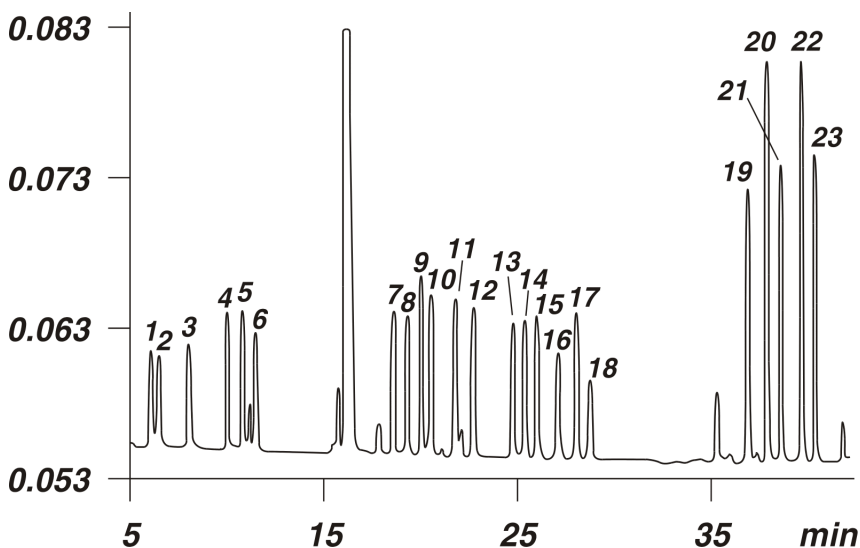

Dabsyl Derivated Phospho Amino Acids

ANALYTES:

- 1** Phosphothreonine
- 2** Phosphoserine
- 3** Phosphotyrosine
- 4** Aspartic acid
- 5** Glutamic acid
- 6** Carboxymethylcysteine
- 7** Serine
- 8** Threonine
- 9** Glycine
- 10** Alanine
- 11** Proline
- 12** Valine
- 13** Methionine
- 14** Isoleucine
- 15** Leucine
- 16** Arginine
- 17** Phenylalanine
- 18** Cysteine
- 19** Ammonia
- 20** Lysine
- 21** Histidine
- 22** Dityrosine
- 23** Tyrosine

