

Methocarbamol and Caffeine (USP XXI 2nd Supplement, p. 1865)

Column: Ultracarb™ 5 µm ODS (30) 60 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-0351-E0

Elution Type: Isocratic

Eluent A: 0.05 M Monobasic potassium phosphate pH 4.6

Eluent B: Methanol

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
1	1	0	70	30

Flow Rate: 1 mL/min


Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 274 nm (ambient)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C18 4 x 3.0mm, 10/Pk Part No.: AJ0-4287

- Holder Part No.: KJ0-4282

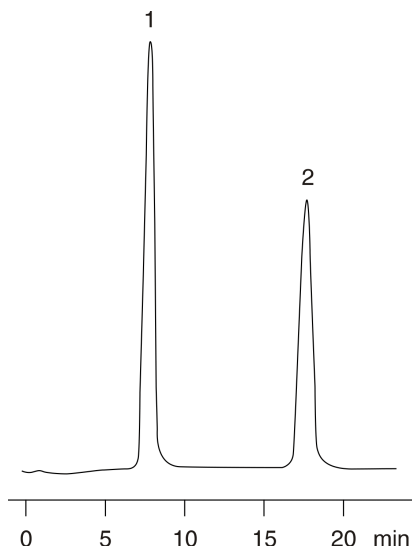

ANALYTES:

1 Caffeine

2 Methocarbamol



Products used in this application:

Ultracarb™

