

Trimethoprim and Sulfamethoxazole

Column: Ultracarb™ 5 µm ODS (30) 60 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-0351-E0

Elution Type: Isocratic

Eluent A: Water

Eluent B: Acetonitrile containing 0.1% TEA, pH 5.9 with acetic acid

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
1	0	80	20	

Flow Rate: 2 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.


- SecurityGuard Cartridges, C18 4 x 3.0mm, 10/Pk Part No.: AJ0-4287

- Holder Part No.: KJ0-4282



Products used in this application:

Ultracarb™

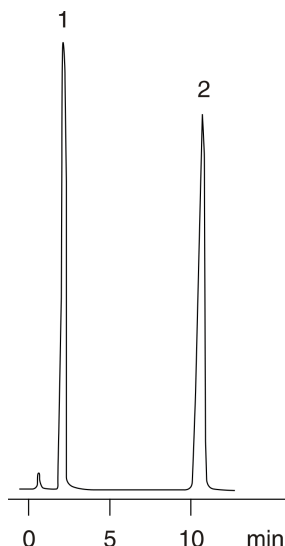

ANALYTES:

- 1 Trimethoprim
- 2 Sulfamethoxazole

