

Wheat Protein Extract - Jupiter vs. Vydac

Column: Jupiter® 5 µm C18 300 Å, LC Column 250 x 4.6 mm, Ea

Dimensions: 250 x 4.6 mm ID

Order No: 00G-4053-E0

Elution Type: Isocratic

Eluent A: Water and acetonitrile (both 0.1% TFA)

Gradient	Step No.	Time (min)	Pct A
Profile:	1	0	100

Flow Rate: 0.5 mL/min


Col. Temp.: ambient

Detection: Unknown @ nm (ambient)


Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Widedpore C18 4 x 3.0mm, 10/Pk Part No.: AJ0-4321

- Holder Part No.: KJ0-4282

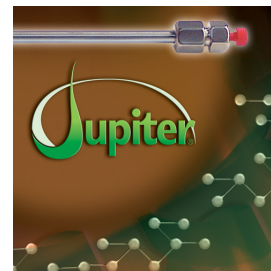

Products used in this application:

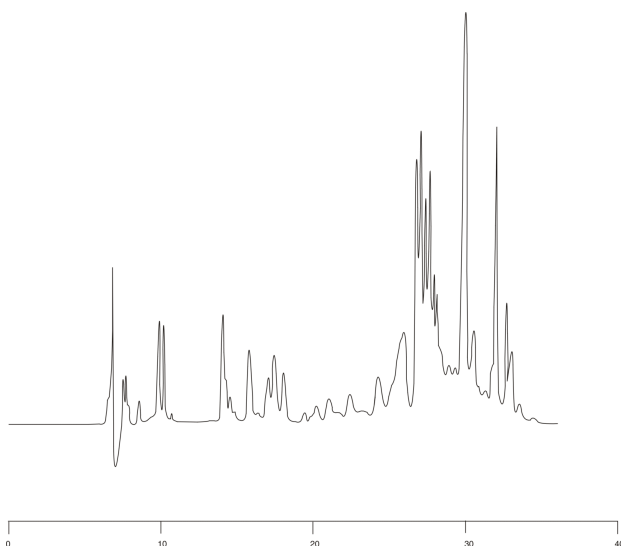

ANALYTES:

1 Wheat protein

