

Tryptic Peptide - Jupiter vs. Vydac

Column: Jupiter® 5 µm C5 300 Å, LC Column 250 x 4.6 mm, Ea

Dimensions: 250 x 4.6 mm ID

Order No: 00G-4052-E0

Elution Type: Isocratic

Eluent A: 0.1% TFA

Eluent B: Acetonitrile/0.1% TFA in water (95:5)

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
1	0	50	50	

Flow Rate: 1 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 214 nm (ambient)


Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Widedpore C5 4 x 3.0mm, 10/Pk Part No.: AJ0-4327

- Holder Part No.: KJ0-4282



Products used in this application:

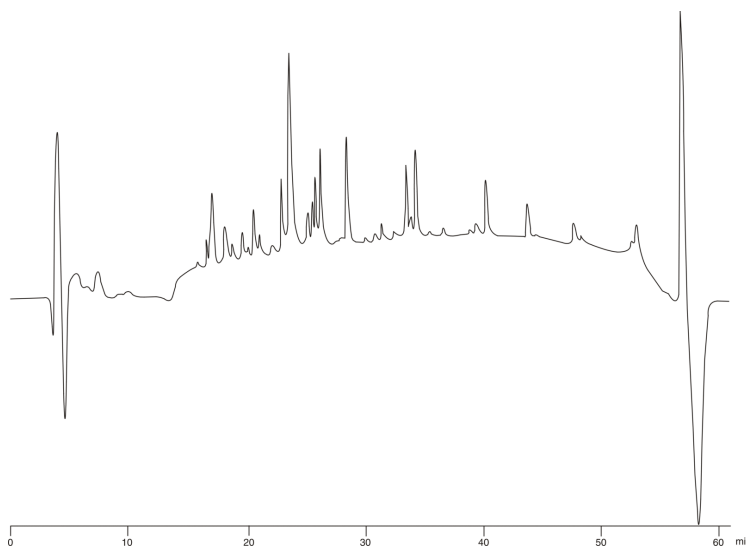

ANALYTES:

1 Protein

