

Doxepin HCl - USP Method


Column: Luna® 5 µm C8(2) 100 Å, LC Column 125 x 4.6 mm, Ea
Dimensions: 125 x 4.6 mm ID
Order No: 00E-4249-E0
Elution Type: Isocratic
Eluent A: 61:39 200mM Monobasic Sodium Phosphate:Methanol, adjusted to pH 2.5
Gradient

Step No.	Time (min)	Pct A
1	0	100

Profile:
Flow Rate: 1 mL/min
Col. Temp.: 50 °C
Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)
Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.
 - SecurityGuard Cartridges, C8 4 x 3.0mm, 10/Pk Part No.: AJ0-4290
 - Holder Part No.: KJ0-4282



Products used in this application:

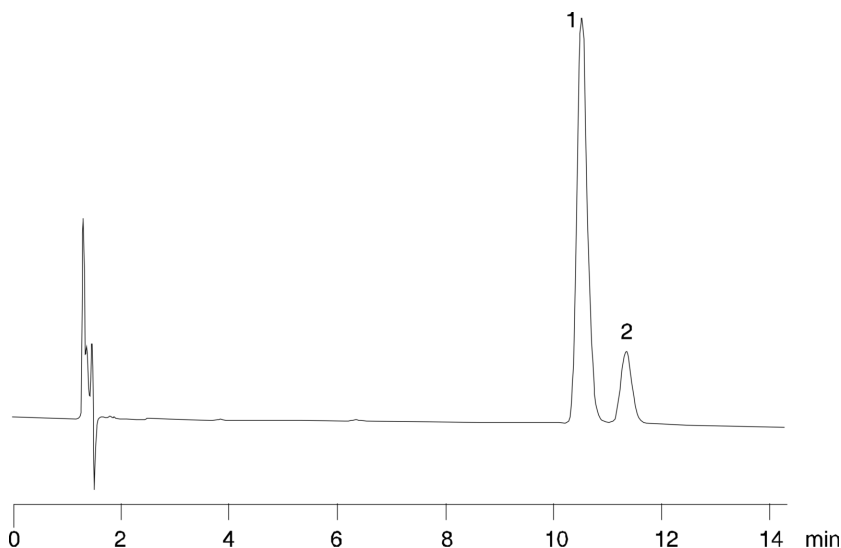

ANALYTES:

- 1 Doxepin
- 2 Doxepin

