

Nortriptyline HCl Capsules - USP Method

Column: Luna® 5 µm CN 100 Å, LC Column 250 x 4.6 mm, Ea

Dimensions: 250 x 4.6 mm ID

Order No: 00G-4255-E0

Elution Type: Isocratic

Eluent A: Acetonitrile

Eluent B: Methanol

Eluent C: 17mM potassium phosphate, pH 6.7

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B	Pct C
1	1	0	40	43	17

Flow Rate: 2.5 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 239 nm (ambient)


Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, CN 4 x 3.0mm, 10/Pk Part No.: AJ0-4305

- Holder Part No.: KJ0-4282



Products used in this application:

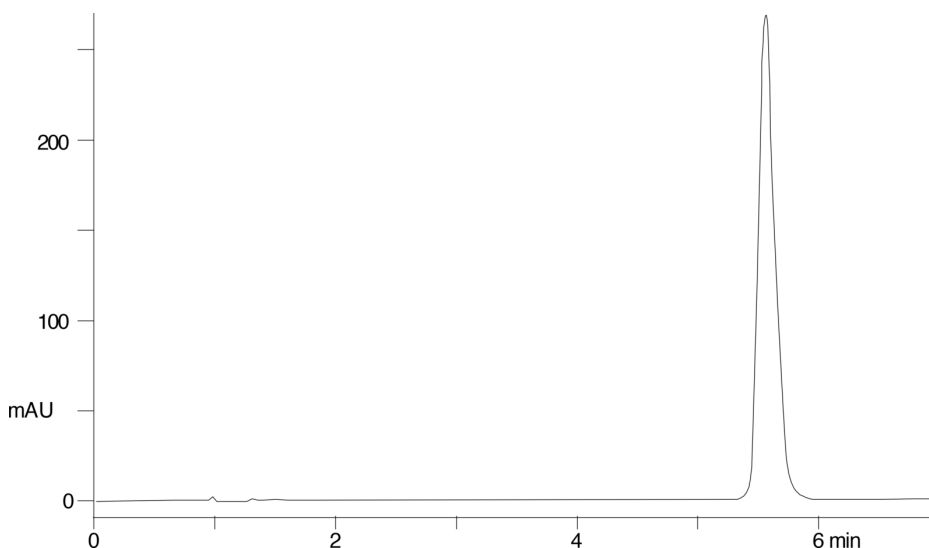

ANALYTES:

1 Nortriptyline

Retention Time: 5.556 min

