

Nortriptyline HCl Capsules - USP Acceptable Modified Method

Column: Luna® 3 µm CN 100 Å, LC Column 100 x 4.6 mm, Ea

Dimensions: 100 x 4.6 mm ID

Order No: 00D-4254-E0

Elution Type: Isocratic

Eluent A: Acetonitrile

Eluent B: Methanol

Eluent C: 17mM phosphate buffer (pH6.7)

Gradient	Step No.	Time (min)	Pct A	Pct B	Pct C
Profile:	1	0	40	43	17

Flow Rate: 1.5 mL/min

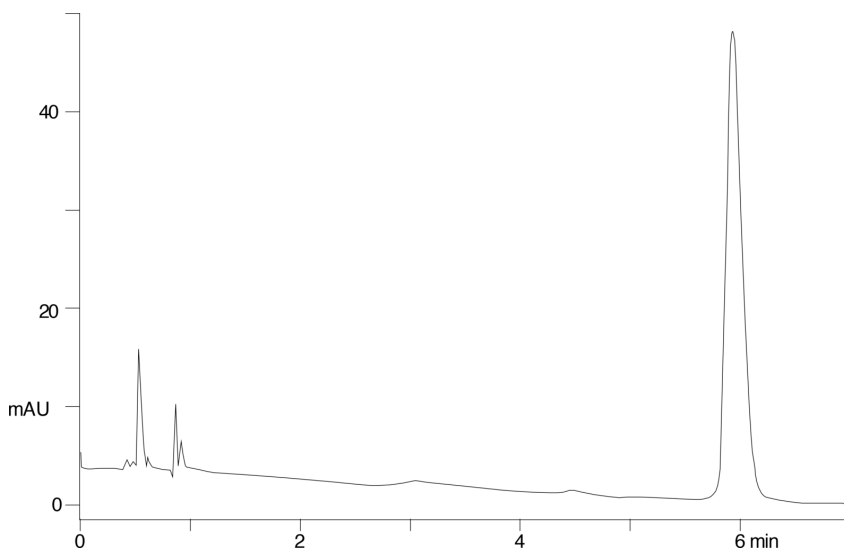
Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 239 nm (ambient)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, CN 4 x 3.0mm, 10/Pk Part No.: AJ0-4305

- Holder Part No.: KJ0-4282

**ANALYTES:**

1 Nortriptyline

Retention Time: 5.928 min



Products used in this application:

