

Quinones

Column: Luna® 3 µm C18(2) 100 Å, LC Column 75 x 4.6 mm, Ea

Dimensions: 75 x 4.6 mm ID

Order No: 00C-4251-E0

Elution Type: Isocratic

Eluent A: CH₃OH

Eluent B: Water

Gradient	Step No.	Time (min)	Pct A	Pct B
----------	----------	------------	-------	-------

Profile:	1	0	89	11
----------	---	---	----	----

Flow Rate: 1 mL/min


Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ nm (ambient)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C18 4 x 3.0mm, 10/Pk Part No.: AJ0-4287

- Holder Part No.: KJ0-4282

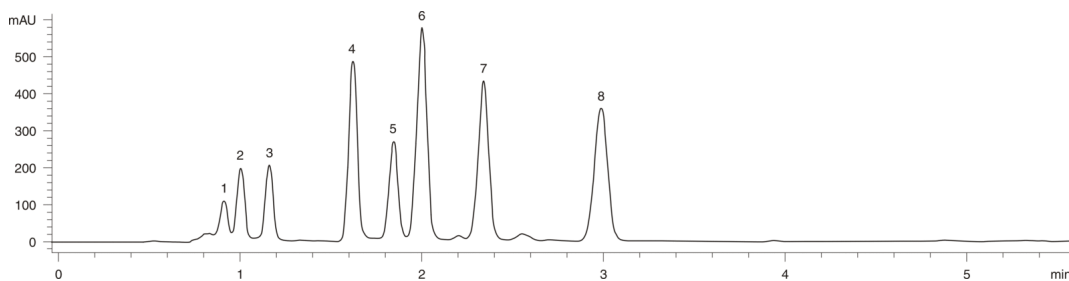

Products used in this application:

**ANALYTES:**

1 Quinones