


Benzalkonium Chloride

Column: Luna® 5 µm CN 100 Å, LC Column 150 x 4.6 mm, Ea
Dimensions: 150 x 4.6 mm ID
Order No: 00F-4255-E0
Elution Type: Isocratic
Eluent A: Acetonitrile
Eluent B: Tetramethyl ammonium chloride buffer 50 mM, pH 3.0
Gradient

Step No.	Time (min)	Pct A	Pct B
1	0	55	45

Profile:
Flow Rate: 1 mL/min
Col. Temp.: ambient
Detection: UV-Vis Abs.-Variable Wave.(UV) @ 214 nm (ambient)
Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.
 - SecurityGuard Cartridges, CN 4 x 3.0mm, 10/Pk Part No.: AJ0-4305
 - Holder Part No.: KJ0-4282

9905

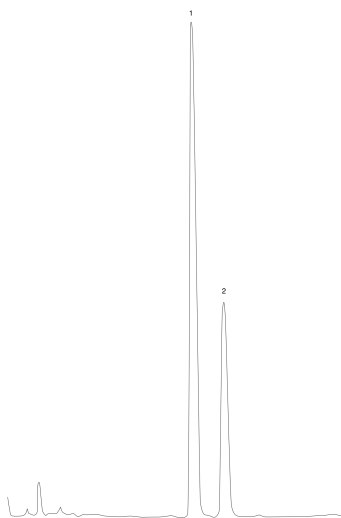

ANALYTES:

- 1 Benzalkonium chloride-C12
- 2 Benzalkonium chloride-C14



Products used in this application:

